



## ***Fractional RF Micro needle Machine***

- ☑ Face lift
- ☑ skin rejuvenation
- ☑ wrinkle
- ☑ scar remover
- ☑ skin tightening
- ☑ acne treatment



# Treatment principle




The patented Flow Needling Technology serially inserts the microneedles row-by-row for a more comfortable procedure




Serial insertion reduces unnatural stretching in order to further reduce pain and downtime. The depth of the coagulation zone is precisely controlled.




Sub-dermal delivery of RF energy eliminates thermal damage to the epidermis providing faster healing time for your patients .

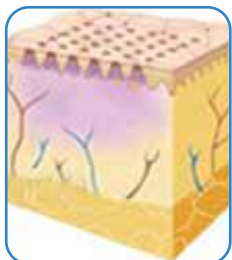


The adjustable depth of the RF energy delivery allows for multiple passes and treatment of delicate areas

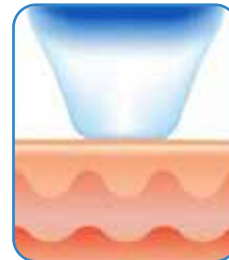


Superficial Fractional RF has two unique channels which combine to provide epidermal ablation and dermal coagulation.

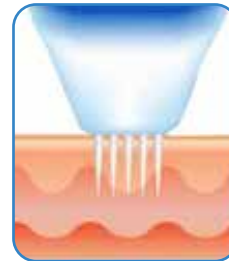

The first channel delivers a deeper zone of controlled thermal effect to activate collagen and skin enhance skin tightening.

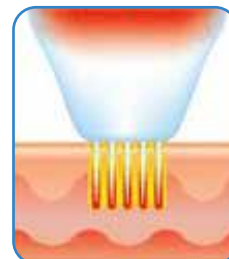

The second channel provides micro-ablation of the upper layers of the epidermis and mild coagulation for controlled fractional ablation and skin rejuvenation with minimally invasive RF delivery.

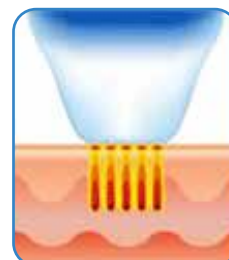

Together , the two channels provide a 3-dimensional treatment zone for superior skin rejuvenation, overall tightening and consistent results you can count on.

# Handle Advantages



MRF-81  
needles



MRF-49  
needles



MRF-25  
needles



SRF

Automatic output controlling ,  
0.3-3mm depth adjustable



# Before & After

## Treatment scope

Skin Tightening, Periorbital Wrinkles, Skin Rejuvenation, Active Inflammatory Acne, Enlarged Pores, Acne Scar, Excessive Sebum Secretion, Hyperhidrosis/Osmidrosis, Stretch Mark



Before



After



Before



After



Before



After



Before



After

# Advantages



- ▶ Ergonomically designed handpiece
- ▶ Precise targeting with adjustable microneedles over 3.0 mm
- ▶ Safely treat all skin types, including tanned skin
- ▶ Specifically target defined depths with precise coagulation
- ▶ Deeper tissue effect than traditional RF and fractional lasers
- ▶ Quick and easy recovery with little downtime

# *Technical parameter*

Frequency	5 MHZ
Total RF Treatment Level	1- 10
Power	80W
Type of needles	81 tips, 49tips, 25tips
Depth of Needle	0.3-3mm(Adjustable)
Switch	Foot Switch
Dimension(CM)	40*38*80
Weight	20kg

## *Contact us*

Beijing Nubway S & T Co., Ltd.

Address: Room 805 No.12 Yabao Road,Chaoyang  
District, Beijing 100020, China

E-mail: [Info@bjnubway.com](mailto:Info@bjnubway.com)  
Phone: 0086-10-51206522

Website: [www.nubway.net](http://www.nubway.net)